N1

User Manual

breo®

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ENGLISH

User Manual

Thank you for purchasing the breo products.
Before using this product,
please read this manual carefully
for convenient use.

Contents

Product Overview and Components	32
Product Features	33
How to Use	34
Product Specifications and Problem Solving	36
Maintenance	37
Cautions	38
Guidance on Product Replacement / Refund	40
Guidance on Quality Assurance	41
Product Warranty	42

Product Overview and Components

Product Overview

N1 will relax your neck through unique acupressure effects that make you feel like you're getting acupressure from an expert. Thanks to its luxurious, knob-detachable design, you can also use the Product for other body parts than neck.

Product Components

N1, User Manual, Charging adapter

N1


User Manual


Charging adapter

Product Features

Product Features

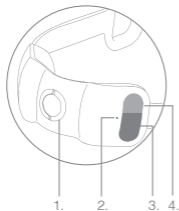1. Width Control Knob
2. LED Indicator
3. Function Control Button (Select Mode button, Power ON/OFF button)

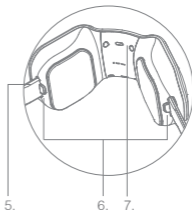4. LCD Display Window
5. Charging Terminal
6. Detach Knob Belt Button
7. IR Lamp

How to Use

Preparations before Use

1. Make sure that the battery has been fully charged.
2. You should always charge the battery with a reliable charging adapter that comes with the Product.

Charging the Product

1. Plug the charging adapter's DC plug into the charging terminal.
2. The LED indicator blinks while charged.
3. The Product will automatically stop while it is being charged.
4. It takes about 4 to 5 hours to charge the Product.

Wearing the Product

Insert the knob belt into N1. For proper treatment, you can adjust the width using the width control knob. Pull up the silver-colored semi-circle (marked "V") in the center of the Product and then turn it left or right for optimal adjustment. Place the relaxer onto the back of your neck and properly position the relaxer while holding the both knob ends.

Operation

If you long-press the Power ON/OFF button, the Product will be activated while the LED indicator blinks in red. The default operating time is 15 minutes, and the default mode is Mode 1.


How to Use

Mode selection

If you short-press the Select Mode button, mode will change in a preset sequence. 4 modes are available.

Mode 1	Mode 2	Mode 3	Mode 4
Left cheek Clockwise rotation ----- Right cheek Counterclockwise rotation	Left cheek Counterclockwise rotation ----- Right cheek Clockwise rotation	Both cheeks Counterclockwise rotation	Both cheeks Clockwise rotation

Optimal effects can be achieved when you regularly use the Product. It is highly recommended to use the Product twice, 15 minutes each.

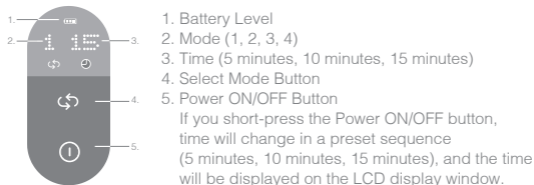Power OFF

If you long-press the Power ON/OFF button, the Product will be deactivated while the LED indicator is turned off. (When the set time elapses, the Product will be automatically turned off.)

Product Specifications and Problem Solving

Product specification

Rated Voltage	DC 5V 1A
Power Consumption	5W
Dimensions / Weight	(W)220mm X (L)105mm X (H)120mm / 835g
Operating Environment	Temperature : 5°C ~ 40°C / Humidity : ≤80% Air pressure : 860hpa ~ 1060hpa

Problem Solving

If you experience any problem using the Product, please check the methods specified below. If the problem persists, turn the Product off and send it to the service center. Never disassemble or repair the Product on

Problems	Possible Causes	Solution
If the Product cannot be turned on	Low or dead battery	Charge the battery.
If rotation is too weak	Low or dead battery	Charge the battery.

Maintenance

Storage precautions

Keep the product in a safe place to avoid exposure to moisture and contamination. If the storage period exceeds 6 months, check the operation of product and turn on the power for about 30 minutes. Avoid the direct sunlight and corrosive gases.

Transporting precautions

This product is designed for seismic resistance and can be transported by plane, train, car, and ship. Avoid snow, rain, and impact as much as possible and put it in the right direction.

Troubleshooting

If a malfunction occurs during operation, check the main battery first. If the trouble cannot be solved, turn off the power and send it to our A/S service center. Do not self-disassemble or repair. If the voltage exceeds the rated range of product, do not turn on the product. Check the cables of power supply periodically.

Cautions

※ Please read the following cautions carefully and accordingly use the Product.

1. Do not use the product after you drink, when you don't feel good, or while you are driving.
2. Any person with a stroke, vascular disease or hypoglycemia is not allowed to use the Product.
3. If a body part to be treated with the Product has a rash or infection, or has undergone a surgery or treatment, you cannot use the Product.
4. Any child aged under 7, any immobile patient, or any physically or mentally challenged person is not allowed to use the Product.
5. Any person with a pacemaker or electronic prosthesis, or any pregnant woman is not allowed to use the Product.
6. The Product does not substitute any diagnosis or treatment by a medical specialist. If you want to confirm that you have no health issue for using the Product, please consult with your doctor.

Cautions

- ※ Please read the following cautions carefully and accordingly use the Product.
7. If a body part to be treated with the Product has a rash or infection, or has undergone a surgery or treatment, it is not recommended to use the Product.
 8. The Product has the heating function. Therefore, any person who is insensitive to heat should take extra care.
 9. If you experience any problem or displeasure while using the Product, stop and turn off the Product immediately, and perform relevant inspection.
 10. In terms of a charging adapter, use the enclosed stabilized one only. Be careful for electrical safety.
 11. The Product is a household product. Therefore, it is strictly prohibited from using the Product for any commercial, treatment or diagnostic purpose.
 12. Use a sterilized cotton dipped in neutral detergent when cleaning the Product. Never use water as cleaning liquid.

Guidance on Product Replacement / Refund

1. If a critical repair due to any performance- or function-related under a normal operating environment is required defect within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have main tained their original conditions before unpacked.

If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center(1588-2336). Batteries are not applicable to any warranty service.

Guidance on Quality Assurance

Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty service is valid for one year from the date of purchase.

Quality Assurance Conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
2. For any defect due to a user's fault(e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty.(e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose(e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster(e.g. typhoon, earthquake, flood, etc.), no warranty service can be provided.

Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

Product Warranty

Name of Product :

Serial Number :

Date Purchased :

Place Purchased :

Name of Customer :

Contact :

Address :

Guidance on Product Replacement/Refund

1. If a critical repair due to any performance- or function-related under a normal operating environment is required defect within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have main tained their original conditions before unpacked.

If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center(1588-2336). Batteries are not applicable to any warranty service.

Product Warranty

Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty service is valid for one year from the date of purchase.

Quality Assurance Conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
2. For any defect due to a user's fault(e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty.(e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose(e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster(e.g. typhoon, earthquake, flood, etc.), no warranty service can be provided.

Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.