H1

User Manual

breo®

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ENGLISH

User Manual

Thank you for purchasing the breo products.
Before using this product,
please read this manual carefully
for convenient use.

Contents

Product Overview and Components	43
Functions	44
How to Use	46
Product Specifications and Troubleshooting	49
Maintenance	50
Cautions	51
Guidance on Product Replacement/Refund	53
Guidance on Quality Assurance	54
Product Warranty	55
App User Manual	57

Product Overview and Components

Product Overview

H1 will alleviate any stress or fatigue by massaging certain points on your palm. Furthermore, it will enhance the blood circulation, hence resulting in healthier organs and body parts.

Products Components

H1, Charging Adapter, User Manual

H1


Charging Adapter

User Manual

Fuctions

Main Body

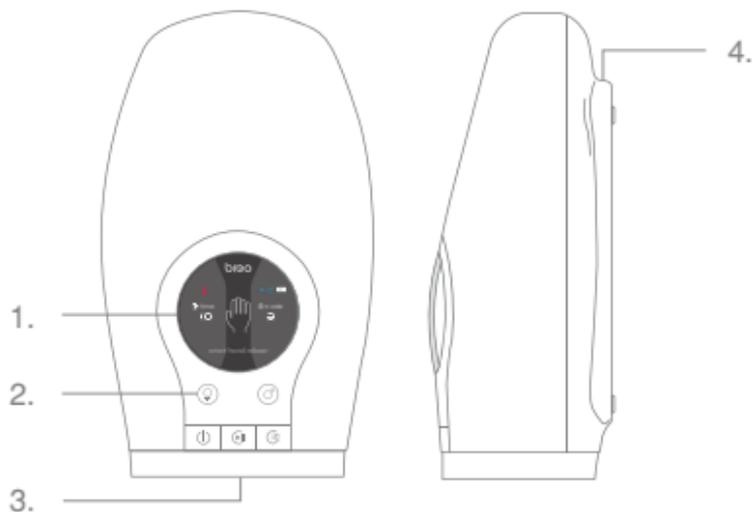1. Display
2. Button
3. Hand Insertion Part
4. Charging Terminal

Fuctions

Display and Button Details

- | | |
|--------------------------------|----------------------|
| 1. Female Mode | 9. Bluetooth |
| 2. Male Mode | 10. Male Mode |
| 3. Power Button | 11. Battery Capacity |
| 4. Start/Pause Button | 12. Method |
| 5. Air Pressure Control Button | 7. Timer |
| 6. Female Mode | 8. Operating Status |

How to Use

Power ON/OFF

1. If you press the Power button, the relaxer will be activated.
2. If you press the Start/Pause button, the relaxer will start running.
3. If you press the Start/Pause button, the relaxer will temporarily stop.
4. If you press the Power button, the relaxer will be deactivated.


Operating the Product

1. Switch to male mode : Press the (♂) button. You can select a mode from Mode 1 to Mode 5.
 2. Switch to female mode : Press the (♀) button. You can select a mode from Mode 1 to Mode 5.
 3. Air pressure control: If air pressure is too strong, press the (🌀) button to depressurize the air.
- ※ For more details, please refer to the “User Manual – Applications” inside the product box.

How to Use

Male Mode

* See the figure below.

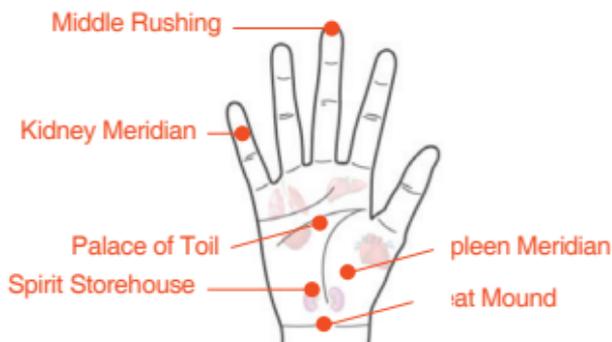Mode1. When you go to the bathroom too frequently

If you suffer from any residual urine or frequent urination, press the "Kidney Meridian" (first knuckle of a little finger) in order to relax your musculoskeletal system and balance your kidney.

Mode2. When you feel stressed out

Whenever you have chest pain or feel stressed out, press the "Palace of Toil" to relax your mind and body.

Mode3. When you are not feeling well

When you feel dizzy or distracted, press the "Spleen Meridian" to enhance the blood circulation in your spleen.

Mode4. When you need a deep sleep Whenever you have a hard time going to bed, there must be some problems associated with the flow of your acupuncture points. Press the "Palace of Toil", "Middle Rushing" and "Great Mound" to get a deep sleep

Mode5. When you feel bloated

Whenever you feel bloated, press the "Palace of Toil" and "Great Mound" to ease your stomach and enhance your blood circulation.

Product Specifications and Troubleshooting

Product Specifications

Rated Voltage	DC 5V	Power Consumption	6W
Bluetooth	Bluetooth 4.0		
Battery Type	Large lithium polymer battery		
Dimensions / Weight	(W)235mmX(L)145mmX(H)98mm / 740g		
Operating Environment	Temperature : 5°C ~ 40°C / Humidity : ≤80% Air pressure : 860hpa ~ 1060hpa		

Troubleshooting

If you experience any problem using the Product, please check the methods specified below. If the problem persists, turn the Product off and send it to the service center. Never disassemble or repair the Product on your own.

Problems	Possible Causes	Solution
If the Product cannot be turned on	Low or dead battery	Charge the battery.
If air pressure or vibration is too weak	Low or dead battery	Charge the battery.
If the Product cannot be charged	Poor charger connection	Make sure the charger has been properly connected.

Maintenance

Storage precautions

Keep the product in a safe place to avoid exposure to moisture and contamination. If the storage period exceeds 6 months, check the operation of product and turn on the power for about 30 minutes. Avoid the direct sunlight and corrosive gases.

Transporting precautions

This product is designed for seismic resistance and can be transported by plane, train, car, and ship. Avoid snow, rain, and impact as much as possible and put it in the right direction.

Troubleshooting

If a malfunction occurs during operation, check the main battery first. If the trouble cannot be solved, turn off the power and send it to our A/S service center. Do not self-disassemble or repair. If the voltage exceeds the rated range of product, do not turn on the product. Check the cables of power supply periodically.

Cautions

※ Please read the following cautions carefully and accordingly use the Product.

1. Read the user manual thoroughly and use the product as instructed in the manual.
2. Do not wear any jewelry (e.g. ring, etc.) when using the product.
3. Do not polish your nail within 1 hour prior to using the product.
4. Turn the product on only after inserting your hands in order to prevent any overpressure-related damage to the air bag.
5. Some highly-sensitive users may feel a slight pain when using the product.
6. It is normal to have some marks on your skin immediately after product. Such marks will usually disappear within 15 minutes.
7. Do not use the product after you drink, when you don't feel good, or while you are driving.
8. Do not use or keep the product in any place where it is humid or hot, or where there is any corrosive gas.
9. Do not keep the product near any electric heater or cooker, or expose it to any sunlight.
10. Any person with a stroke, vascular disease or hypoglycemia is not allowed to use the product.
11. Any child aged under 7, any immobile patient, or any physically or mentally challenged person is not allowed to use the product.

Cautions

12. Any person with a pacemaker or electronic prosthesis, or any pregnant woman is not allowed to use the product.
13. If a body part to be treated with the product has a rash or infection, or has undergone a surgery or treatment, you cannot use the Product.
14. The product does not substitute any diagnosis or treatment by a medical specialist. If you want to confirm that you have no health issue for using the Product, please consult with your doctor.
15. If you experience any problem or displeasure while using the product, stop and turn off the product immediately, and perform relevant inspection.
16. You should pull your hands after turning the Product off. If the product suddenly stops running, try pulling your hands by twisting them left and right.
17. The Product is a household product. Therefore, it is strictly prohibited from using the product for any commercial, treatment or diagnostic purpose.
18. Use a sterilized cotton dipped in neutral detergent when cleaning the product. Never use water as cleaning liquid.

Guidance on Product Replacement / Refund

1. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have maintained their original conditions before unpacked.

If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center (1588-2336). Batteries are not applicable to any warranty service.

Guidance on Quality Assurance

Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty period is one year from the date of purchase.

Warranty conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty (e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster (e.g. typhoon, flood, earthquake, etc.), no warranty service can be provided.

※ Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

Product Warranty

Name of Product : _____

Serial Number : _____

Date Purchased : _____

Place Purchased : _____

Name of Customer : _____

Contact : _____

Address : _____

Guidance on Product Replacement/Refund

1. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have maintained their original conditions before unpacked. If you intend to make such request or have any related inquiry, please contact your buyer or our customer center (1588-2336).

Product Warranty

Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty period is one year from the date of purchase.


Warranty conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
 2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
 3. Additional expense may incur for any supplementary services not specified in the product warranty (e.g. delivery service, etc.)
 4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
 5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
 6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
 7. For any defect due to a natural disaster (e.g. typhoon, flood, earthquake, etc.), no warranty service can be provided.
- ※ Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

How to Use

Female Mode

* See the figure below.

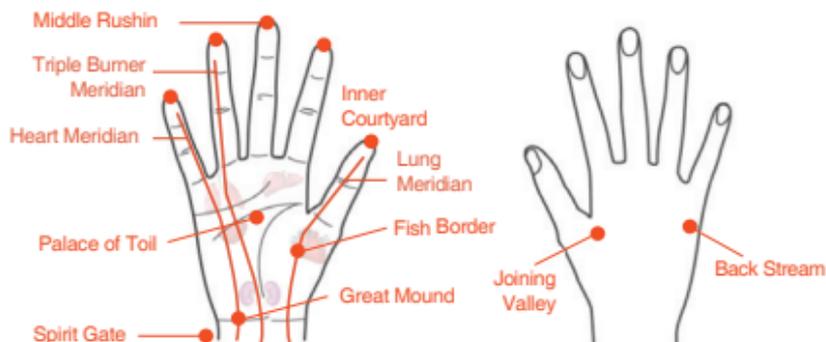Mode 1. When your hands get too dry

For your hands exposed to the external environment (e.g. dust, frequent washing, etc.), massage the "Lung Meridian", "Heart Meridian" and "Triple Burner Meridian" to relax the muscle and skin in your hands.

Mode 2. When you feel cold in your hands or feet

When you feel cold in your hands or feet, or look pale, massage the "Inner Courtyard" (fingertip).

Mode 3. When you have a dark skin tone When you look tired or have dark circles around your eyes, massage the "Palace of Toil", "Fish Border" and "Spirit Gate" to moisten your skin.

Mode 4. When you need a deep sleep Whenever you have a hard time going to bed, there must be some problems associated with the flow of your acupuncture points. Press the "Palace of Toil", "Middle Rushing" and "Great Mound" to get a deep sleep.

Mode 5. When you feel too tight in your shoulder muscle


When you feel too tight in your neck after long hours of smart-phone and computer uses, massage the "Joining Valley" and "Back Stream" to relax your shoulder muscle.

How to Use APP

1. Download ihand application.
(from Google Playstore on Android or App Store for iOS)
2. Turn on iHand relaxer.


3. Launch iHand application on your smartphone.
4. Approve Bluetooth activation
(if Bluetooth is turned on, skip to the next step)

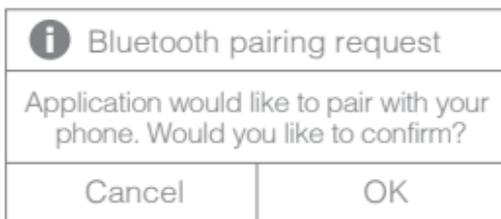
How to Use APP

5. Select iHand from searched devices for connection.
If your smartphone cannot find the device, please check to see if iHand is turned on.

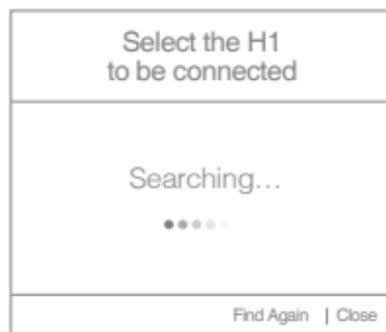
6. Register iHand on Bluetooth.
(When registered on Bluetooth, steps 4-6 may be skipped for further uses.)

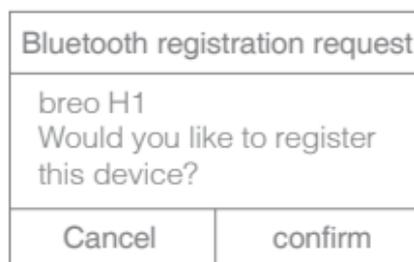
How to Use APP

7. Hold the device in hand, and select mode.

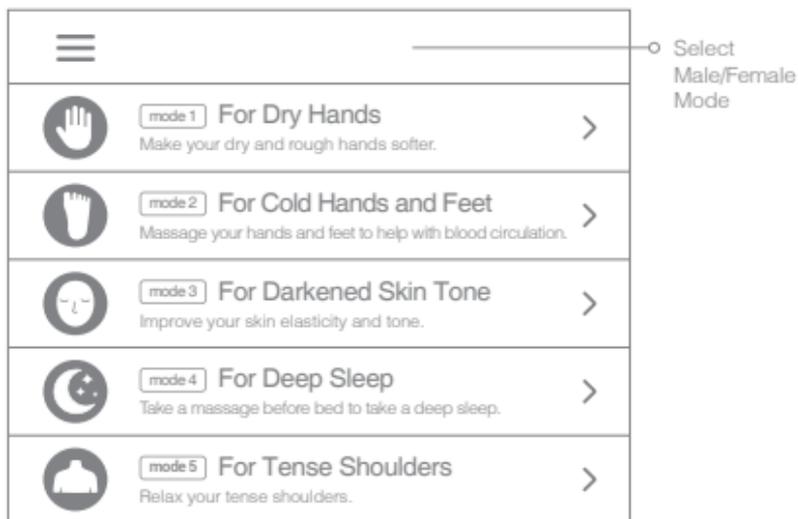8. Press “On” button to begin the massage (when the pressure is too high or too low, change the strength).

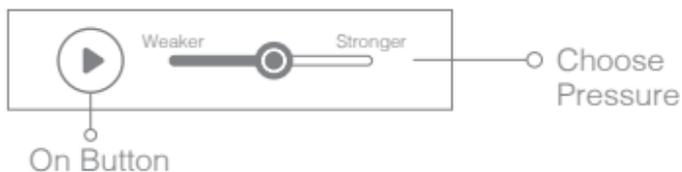9. When done with massage, press “Stop” button to stop the device.