# B1

User Manual

breo®

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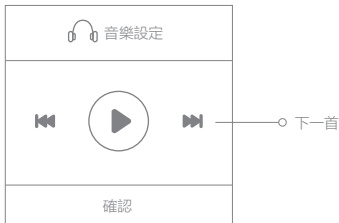
# ENGLISH

User Manual

Thank you for purchasing the breo products.  
Before using this product,  
please read this manual carefully  
for convenient use.

## Contents

---

Product Overview and Components	36
Functions	37
How to Use	38
Product Specifications and Troubleshooting	41
Maintenance	42
Cautions	43
Guidance on Product Replacement/Refund	45
Guidance on Quality Assurance	46
Product Warranty	47
APP User Manual	

## Product Overview and Components

---

### Product Overview

B1 is the world's first smart head&eye relaxer. It massages certain acupuncture points in the head and area around the eyes through air pressure control, acupressure, vibration and heat. In addition, it keeps users healthy by supplying oxygen, improving blood circulation, easing stress and bringing good sleep. In addition, it is designed to adjust a head size conveniently for diverse people. Furthermore, the part around the eyes is detachable so that the Product is available in diverse environments. It presents smart design. You are able to download an app by connecting the Product to your smart phone. In addition, the functions are easy to use so that you can enjoy the massage in comfort.

### Product Components

B1, Charging adapter, Earplug, User manual

B1

Charging adapter

Earplug


User manual

## Fuctions

---

### Fuctions

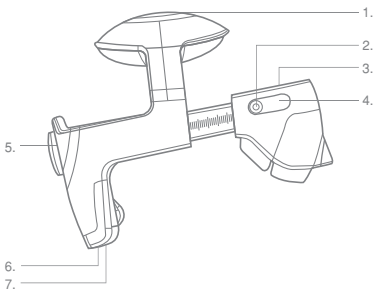1. Height Adjustment  
Button

2. Power Button  
Mode Button

3. Eye Separation Button

4. State Indicator

5. Head Size Adjustment Wheel

6. Earplug Port

7. Charging Port (USB Port)

## How to Use

---

### Preparations before Use

Check the Product before use as follows:

1. Check the device and its components.
2. Check if the battery is sufficiently charged.
3. Check if the charging adapter is free of defect.

### Charge

Connect the Charging Adapter

1. Open the charge cover on the back of the Product.
2. Plug the charging adapter into the power port.
3. If the charge level indicator icon blinks, the battery is being charged. If the bar is fully filled, normal charge is completed.

#### CAUTION 1

Make sure that the Product is turned off during charge.

The normal charge time is 4 hours.

#### CAUTION 2

Do not use the Product during charge.

(The Product does not work during charge. If attempted to charge during the use of the Product, the auto-protection program is initiated and stops the system automatically.)

### Put On

1. Adjust the head size by turning the head size adjustment wheel before putting on the Product.
2. Adjust the height by pressing and holding the height adjustment button short or long.

## How to Use

---

### Switch On

If the power button is pressed and hold short while the Product is off, it is turned on with a beep sound. Basically, 'air pressure + heating + music + 15 min.' mode is set.

### Mode Selection

	MODE1	MODE2	MODE3	MODE4	MODE5
Air Pressure	O	O	O	O	-
Heating	O	O	-	-	-
Time(Min)	15	15	15	15	15
Scalp Massage	-	-	-	-	O
Music	O	O	O	O	O

### Play Music

Once the Product is switched on, music starts.


### Height Adjustment

Press and hold the height adjustment button short to raise the height, and press and hold the button long to lower the height. (This button is active in all modes.)

## How to Use

---

\* See the figure below.

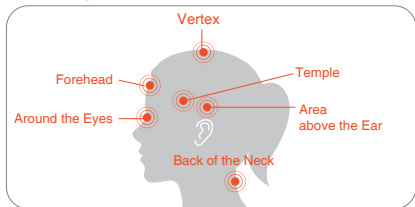**Mode 1. If stressed** If you feel not good because of stress, massage your head, neck and area around the eyes, promote blood circulation and ease tense muscles.

**Mode 2. Tired eyes** If you feel dry and uncomfortable in the eyes, massage the area around the eyes and part above your ears, which would in turn promote blood circulation, provide oxygen and relive fatigue in the eyes.

**Mode 3. Stiff neck** If you feel stiff in your neck because a long use of your smartphone or desktop computer, massage your forehead and back of your neck and relax tense muscles, refreshing your head.

**Mode 4. Headache** If physical and mental stresses are accumulated, you could suffer from a headache. Massage the temple and area around the eyes to relieve the fatigue and ease pain.

**Mode 5. I need concentration** If you feel less energy and find it hard to concentrate, massage the top of your head. It promotes blood circulation and supplies oxygen to the brain, improving concentration and memory.

## Product Specifications and Troubleshooting

### Product Specifications

Rated Voltage	DC 5V
Power Consumption	5W
Dimensions / Weight	(W)210mmX(L)254mmX(H)242mm / 1075g
Operating Environment	Temperature : 5°C ~ 40°C / Humidity : ≤80% Air pressure : 860hpa ~ 1060hpa

### Troubleshooting

If you experience any problem using the Product, please check the methods specified below. If the problem persists, turn the Product off and send it to the service center. Never disassemble or repair the Product on your own.

Problems	Possible Causes	Solution
If the Product cannot be turned on	Low or dead battery	Charge the battery.
If air pressure or vibration is too weak	Low or dead battery	Charge the battery.
If the Product cannot be charged	Poor charger connection	Make sure the charger has been properly connected.

## Maintenance

---

### Storage precautions

Keep the product in a safe place to avoid exposure to moisture and contamination. If the storage period exceeds 6 months, check the operation of product and turn on the power for about 30 minutes. Avoid the direct sunlight and corrosive gases.

### Transporting precautions

This product is designed for seismic resistance and can be transported by plane, train, car, and ship. Avoid snow, rain, and impact as much as possible and put it in the right direction.

### Troubleshooting

If a malfunction occurs during operation, check the main battery first. If the trouble cannot be solved, turn off the power and send it to our A/S service center. Do not self-disassemble or repair. If the voltage exceeds the rated range of product, do not turn on the product. Check the cables of power supply periodically.

## Cautions

---

※ Please read the following cautions carefully and accordingly use the Product.

1. Read the User Manual carefully and follow the instructions.
2. Do not use the Product if you had alcohol, feel uncomfortable or are driving a car.
3. Do not use, keep or charge the Product in a wet or hot place or the area with corrosive gas around.
4. Do not leave the Product around the electric heater or cooking utensil. Avoid a direct sunlight.
5. The Product is available for the adults with 56-64cm in head circumference only.
6. If you suffer from palsy, stroke, blood pressure disease, hypoglycemia, eye disease(e.g., glaucoma, cataract, detached retina, etc.)and brain disease, do not use the Product.
7. For young children (7 years old or under), patients with limited mobility and physically/mentally disabled people, use the Product in the presence of a caregiver. If any discomfort is felt during the use of the Product, stop the Product immediately.
8. If you carry any medical device such as pacemaker and other medical aids on their body or pregnant, do not use the Product.
9. Take off your hair accessories, contact lens and glasses before the use of the Product.
10. If you had an eye surgery such as LASIK, LASEK or lens implant surgery, do not use the Product. If you insist to use it, make sure that you consult with your doctor first.

## Cautions

---

11. The Product does not substitute a professional's diagnosis or treatment. If necessary, consult with a professional before using the Product.
12. If there is any abnormal symptom such as rash or inflammation, or you had a surgical operation and are in the middle of treatment in the target area, do not use the Product.
13. The heating function has nothing to do with hair loss, atopic dermatitis, acne and allergy.
14. The mode setting becomes active once you put on the Product. Otherwise, it could be damaged by exceeding the air pressure limit.
15. If you feel uncomfortable during use, stop the Product immediately and cut off the power. Then, check the Product.
16. For a charging adaptor, use the enclosed DC5V adapter only. Be careful for electrical safety.
17. Turn on the power after putting on the Product. Do not lean on a wall or walk with the Product on.
18. The Product is designed for home use only. Do not use it for commercial or any other (e.g., medical treatment or diagnosis, etc.) purposes.
19. For cleaning, use the sterilized cotton ball with a neutral detergent. Do not EVER try to clean the Product with water.
20. When the heating function is used, if you are dull in feeling temperature, be careful during use.

## Guidance on Product Replacement / Refund

---

1. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have maintained their original conditions before unpacked.

If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center (1588-2336). Batteries are not applicable to any warranty service.

## Guidance on Quality Assurance

---

### Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty period is one year from the date of purchase.

### Warranty conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty (e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster (e.g. typhoon, flood, earthquake, etc.), no warranty service can be provided.

※ Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

## Product Warranty

---

Name of Product : \_\_\_\_\_

Serial Number : \_\_\_\_\_

Date Purchased : \_\_\_\_\_

Place Purchased : \_\_\_\_\_

Name of Customer : \_\_\_\_\_

Contact : \_\_\_\_\_

Address : \_\_\_\_\_

### Guidance on Product Replacement/Refund

1. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function-related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have maintained their original conditions before unpacked. If you intend to make such request or have any related inquiry, please contact your buyer or our customer center (1588-2336).

## Product Warranty

---

### Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty period is one year from the date of purchase.

### Warranty conditions



1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
  2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
  3. Additional expense may incur for any supplementary services not specified in the product warranty (e.g. delivery service, etc.)
  4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
  5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
  6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
  7. For any defect due to a natural disaster (e.g. typhoon, flood, earthquake, etc.), no warranty service can be provided.
- ※ Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

## How to Use

---

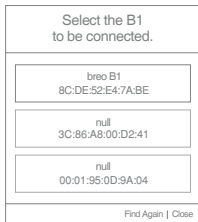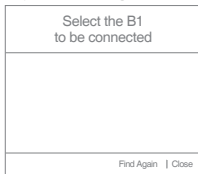1. Please download B1 APP  
( through Android play store / iOS App Store)
2. Put on B1 and start the massage by turning on the power.
3. Open the bluetooth in the mobile phone set and select the "B1" in the bluetooth device list searched.
4. Connect the B1 to the bluetooth of the mobile phone.
5. When the product is successfully connected to the bluetooth, the builtin music will stop. (If the built-in music is still in play, please check if the bluetooth of the mobile phone is successfully connected)
6. Run the B1 APP installed on the phone.
7. When the connection choosing image appears, connect the B1 again in the bluetooth device list searched. If the clues of "can not find the equipment" appears:
  - 1) Please confirm whether the phone bluetooth is connected with the B1
  - 2) Click "find again"

## How to Use APP

---

- If the APP bluetooth is not successfully connected, click on the bluetooth button in the APP and select the B1 in the bluetooth list appeared.
- After disconnecting the bluetooth, it takes 5 seconds to reconnect.
- The B1 is successfully connected to the phone bluetooth.(after connecting the phone bluetooth, the steps of 3-5 could be ignored when reusing the product)
- Choose the required massage mode.

massage mode	
 mode 1	If stressed →Message your head and neck and ease tense muscles.
 mode 2	Tired eyes →Get rid of the fatigue in the eyes through massage.
 mode 3	Stiff neck →Soothe and relax tense muscles in the neck.
 mode 4	Headache →Massage the temple and ease the pain.
 mode 5	I need concentration →Clear your head through massage.

## How to Use APP


---

12. Adjust the massage strength.

13. The bluetooth connection and music playback can not be executed under battery symbol flashing state and low electricity condition.


14. Press  and set the timer.
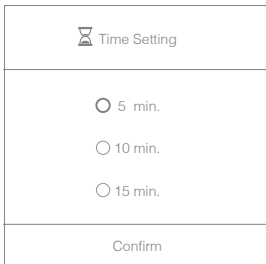
## How to Use APP

---

15. Press  and set the heating function.

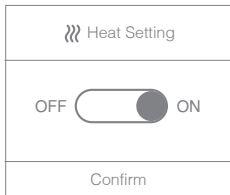16. Press  and set the music.