F1

User Manual

breo®

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ENGLISH

User Manual

Thank you for purchasing the breo products.
Before using this product,
please read this manual carefully
for convenient use.

Contents

Product introduction and configuration	25
Configuration of product functions	26
How to use	27
Product specifications	27
Troubleshooting	28
Maintenance	28
Precautions	29
Information for product exchange / refund	30
Product warranty information	31
Product warranty	32

Product introduction and configuration

Product introduction

F1 loosens the tired foot with the acupressure that allows people to feel the touch of a person

Product configuration

F1, User manual and Charging adapter

F1


User Manual

Charging Adapter

Configuration of product functions

Configuration of product functions

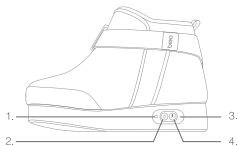- | | |
|--------------------------|-----------------|
| 1. DC charging terminal | 3. Status light |
| 2. Mode selection button | 4. Power button |

Product basic functions

Air pressure	Time
Strong and weak cycling	15 minutes

How to use

Preparations before use

1. Make sure that the battery is fully charged.
2. After wearing the product, turn on the power and start operation.
This prevents the product from being damaged by excessive air pressure when it is operated first before the wearing.

Operation

Power on: When pressing the power button, operation is started after the "beep" sound.

Power off: When pressing the power button, operation is terminated after the "beep, beep, beep" sounds.

How to use

Mode

Mode selection: At the start of operation, it goes to mode 1, and when the mode selection button is pressed, it switches to the next mode after the beep.

Mode type: When mode is changed, status light blinks according to the mode.

- Mode 1: Air pressure (weak) + air acupressure (strong) for 15 minutes / cycle
Status light flashes red one time per second.
- Mode 2: Air pressure (weak) for 15 minutes
Status light flashes red two times every two seconds.
- Mode 3: Air pressure (strong) for 15 minutes
Status light flashes red three times every three seconds.

Charging

Charging starts: The beep sounds twice and charging starts.

Charging in progress: Status light blinks red.

Charging complete: The beep sounds twice and status light turns green. (Takes 6-8 hours)

Mode switching

When operating, mode 1 is activated (15 minutes), and the mode is switched by pressing the mode button briefly.

Product specifications

Product specifications

Rated voltage	DC 5V 2A	Power consumption	10W
Battery type	3.7V Lithium battery		
Size / Weight	(W)310mmX(L)125mmX(H)200mm / 1150g ± 30g		
Using environment	Temperature: 5°C ~ 40°C / Humidity: ≤80% Atmospheric pressure: 860hpa ~ 1060hpa		

Troubleshooting

Troubleshooting

If any problems are encountered while using the product, please check with the methods below. If neither of these methods work, turn off the power and leave it to the service center. Please note that if a non-specialist disassembles or repairs the product at will, it will not be possible to provide further services.

Problem	Possible cause	Solution
When the power is not turned on	Insufficient charge or discharge	Please charge.
When not to be charged	Poor charger connection	Make sure that the charger is properly connected.
When the air pressure is weakened	Insufficient charge or discharge	Please charge.

Maintenance

Storage precautions

Keep the product in a safe place to avoid exposure to moisture and contamination. If the storage period exceeds 6 months, check the operation of product and turn on the power for about 30 minutes. Avoid the direct sunlight and corrosive gases.

Transporting precautions

This product is designed for seismic resistance and can be transported by plane, train, car, and ship. Avoid snow, rain, and impact as much as possible and put it in the right direction.

Troubleshooting

If a malfunction occurs during operation, check the main battery first. If the trouble cannot be solved, turn off the power and send it to our A/S service center. Do not self-disassemble or repair. If the voltage exceeds the rated range of product, do not turn on the product. Check the cables of power supply periodically.

Precautions

*** Please read the user's precautions carefully and use the product in a safe and correct manner.**

1. Do not use if you have experienced foot surgery.
2. Do not use if you have the diseases such as stroke, hypertension, hypotension, etc.
3. This product cannot be used by the children under the age of 7, patients with mobility problems, or persons with physical or mental disabilities.
4. This product cannot be used by the persons who have artificial heart pacemakers or various electronic aids and pregnant women.
5. Do not use the product after drinking alcohol, when not feeling well, or while driving.
6. If you weigh more than 100 kilos, do not walk after wearing the product.
7. If you have a foot disease we do not recommend using the product until you are diagnosed in the stabilization period by the healthcare professional.
8. The diagnosis and treatment of medical professionals cannot be replaced by this product. Before using the product, if it is necessary to check whether there is no problem using the product, be sure to consult with the doctor.
9. Keep this product out of the reach of children.
10. If the product is unusual or uncomfortable during use, discontinue the use and turn off the power for checking.
11. Please use only the adapter that comes with the product as the charging adapter and be careful about electric safety.
12. Put on the product, turn on the power, and start using.
13. As this product has been developed as a household product, it shall not be used for other commercial uses, treatments or diagnosis.
14. When cleaning the product, wipe it with the sanitized cotton wool moistened with neutral detergent. Be careful not to let the product come into contact with water.

Information for product exchange / refund

1. If you need important repairs due to the defects in performance or function in normal use within 10 days from the date of purchase, you can receive a product replacement or refund using your purchase receipt.
2. If you need important repairs due to the defects in performance or function in normal use within one month from the date of purchase, you can exchange or repair the product. However, the product shall not be damaged due to other factors other than the product use defects, and can be exchanged, refunded, or repaired only if it is the same as the original shipment including the user manual and other accessories.

If you would like to receive the services mentioned above, or if you have any questions, please contact the place of purchase or customer center (1588-2336).

Product warranty information

Warranty information

The warranty gives you easy access to the product maintenance services. There is only one warranty for each product, and please note that it cannot be reissued when lost.

※The warranty period is one year from the date of purchase.

Warranty conditions

1. You can receive the repair service with the product warranty for one year from the date of purchase.
If the serial number on the product becomes unclear or the serial number and the information registered on the warranty are different, it may be difficult to provide the appropriate service. (The battery is not included in the warranty.)
 2. If errors are occurred due to the user's negligence such as wrong voltage use or inadvertent accident or due to the repair of this product from an unauthorized company, these cannot be guaranteed.
 3. For the additional services or transportation services other than the specified warranty services, costs may be incurred.
 4. For the problems caused by the mistake or negligence of user, problems caused by the disassembly of non-specialists, faults occurred without warranty, and faults caused by force majeure, the quality assurance service will not be provided.
 5. Please enclose the product warranty when applying for the quality assurance service. We are responsible for the collection of used products, and the removal, disposal and recycling of discarded products.
 6. Warranty is not available when used for the purposes other than household (commercial and public facilities, etc.).
 7. Warranty is not available in the event of the natural disasters such as typhoon, earthquake or flood.
- ※ Caution: Please be advised that the details can be modified without notice as the service is continually upgraded.

Product warranty

Product name: _____

Serial number: _____

Purchase date: _____

Point of purchase: _____

Customer name: _____

Contact information: _____

Address: _____

Information related with product exchange / refund

1. If you need important repairs due to the defects in performance or function in normal use within 10 days from the date of purchase, you can receive a product replacement or refund using your purchase receipt.
2. If you need important repairs due to the defects in performance or function in normal use within one month from the date of purchase, you can exchange or repair the product. However, the product shall not be damaged due to other factors other than the product use defects, and can be exchanged, refunded, or repaired only if it is the same as the original shipment including the user manual and other accessories.

If you would like to receive the services mentioned above, or if you have any questions, please contact the place of purchase or customer center (1588-2336). The battery is not included in the warranty.

Warranty information

The warranty gives you easy access to the product maintenance services. There is only one warranty for each product, and please note that it cannot be reissued when lost.

※ The warranty period is one year from the date of purchase.

Warranty conditions

1. You can receive the repair service with the product warranty for one year from the date of purchase. If the serial number on the product becomes unclear or the serial number and the information registered on the warranty are different, it may be difficult to provide the appropriate service. (The battery is not included in the warranty.)
 2. If errors are occurred due to the user's negligence such as wrong voltage use or inadvertent accident or due to the repair of this product from an unauthorized company, these cannot be guaranteed.
 3. For the additional services or transportation services other than the specified warranty services, costs may be incurred.
 4. For the problems caused by the mistake or negligence of user, problems caused by the disassembly of non-specialists, faults occurred without warranty, and faults caused by force majeure, the quality assurance service will not be provided.
 5. Please enclose the product warranty when applying for the quality assurance service. We are responsible for the collection of used products, and the removal, disposal and recycling of discarded products.
 6. Warranty is not available when used for the purposes other than household (commercial and public facilities, etc.).
 7. Warranty is not available in the event of the natural disasters such as typhoon, earthquake or flood.
- ※ Caution: Please be advised that the details can be modified without notice as the service is continually upgraded.