iBrain

User Manual

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ENGLISH

User Manual

Thank you for buying a **breo** product.
Please read this user manual carefully before using
the Product, and keep it in a proper place for future reference.

* You can download Chinese(Simplified, Traditional) and
Japanese User manual from the breo website.

Table of Contents

Product and Components	3
Functions	3
How to Use	4
Cautions	5
Problem Solving	6
Maintenance	7
Guidance on Quality Assurance	8
Product Warranty	9

Product and Components

About the Product

iBrain relaxer is the world's first smart head&eye relaxer. It massages certain acupuncture points in the head and area around the eyes through air pressure control, acupressure, vibration and heat. In addition, it keeps users healthy by supplying oxygen, improving blood circulation, easing stress and bringing good sleep. In addition, it is designed to adjust a head size conveniently for diverse people. Furthermore, the part around the eyes is detachable so that the Product is available in diverse environments. It presents smart design. You are able to download an app by connecting the Product to your smart phone. In addition, the functions are easy to use so that you can enjoy the massage in comfort.

Components

iBrain relaxer, charging adapter (USB port), earplug, user manual (warranty certificate included)

iBrain relaxer

charging adapter(USB port)


earplug

user manual

Functions

Functions

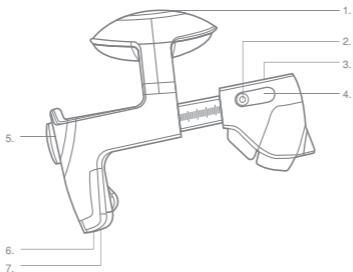- | | | |
|-----------------------------|-------------------------------|-----------------------------|
| 1. Height Adjustment Button | 4. State Indicator | 7. Charging Port (USB Port) |
| 2. Power Button Mode Button | 5. Head Size Adjustment Wheel | |
| 3. Eye Separation Button | 6. Earplug Port | |

How to Use (1)

Preparations Before Use

Check the Product before use as follows:

1. Check the device and its components.
2. Check if the battery is sufficiently charged.
3. Check if the charging adapter is free of defect.

Charge

Connect the Charging Adapter

1. Open the charge cover on the back of the Product.
2. Plug the charging adapter into the power port.
3. If the charge level indicator icon blinks, the battery is being charged. If the bar is fully filled, normal charge is completed.

CAUTION 1. Make sure that the Product is turned off during charge. The normal charge time is 4 hours.

CAUTION 2. Do not use the Product during charge. (The Product does not work during charge. If attempted to charge during the use of the Product, the auto-protection program is initiated and stops the system automatically.)

Put On

1. Adjust the head size by turning the head size adjustment wheel before putting on the Product.
2. Adjust the height by pressing and holding the height adjustment button short or long.

Switch On

If the power button is pressed and hold short while the Product is off, it is turned on with a beep sound. Basically, 'air pressure + heating + music + 15 min.' mode is set.

Mode Selection

	mode1	mode2	mode3	mode4	mode5
Air Pressure	O	O	O	O	-
Heating	O	O	-	-	-
Time(Min)	15	15	15	15	15
Scalp Massage	-	-	-	-	O
Music	O	O	O	O	O

Play Music


Once the Product is switched on, music starts.

Height Adjustment

Press and hold the height adjustment button short to raise the height, and press and hold the button long to lower the height. (This button is active in all modes.)

How to Use (2)

* See the figure below.

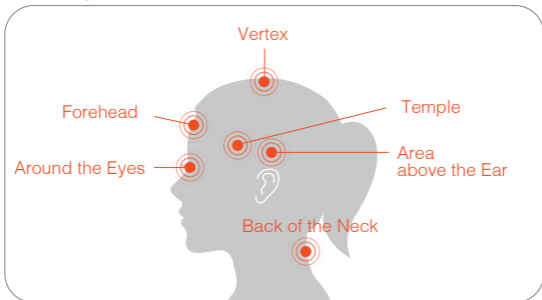Mode 1. If stressed

If you feel not good because of stress, massage your head, neck and area around the eyes, promote blood circulation and ease tense muscles.

Mode 2. Tired eyes

If you feel dry and uncomfortable in the eyes, massage the area around the eyes and part above your ears, which would in turn promote blood circulation, provide oxygen and relive fatigue in the eyes.

Mode 3. Stiff neck

If you feel stiff in your neck because a long use of your smartphone or desktop computer, massage your forehead and back of your neck and relax tense muscles, refreshing your head.

Mode 4. Headache

If physical and mental stresses are accumulated, you could suffer from a headache. Massage the temple and area around the eyes to relieve the fatigue and ease pain.

Mode 5. I need concentration

If you feel less energy and find it hard to concentrate, massage the top of your head. It promotes blood circulation and supplies oxygen to the brain, improving concentration and memory.

Cautions (1)

* Read the CAUTION carefully for the safe and correct use of the Product.

1. Do not use the Product if you had alcohol, feel uncomfortable or are driving a car.
2. The Product is available for the adults with 56-64cm in head circumference only.
3. If you suffer from palsy, stroke, blood pressure disease, hypoglycemia, eye disease (e.g., glaucoma, cataract, detached retina, etc.) and brain disease, do not use the Product.
4. For young children (7 years old or under), patients with limited mobility and physically/mentally disabled people, use the Product in the presence of a caregiver. If any discomfort is felt during the use of the Product, stop the Product immediately.

Cautions (2)

* Read the CAUTION carefully for the safe and correct use of the Product.

5. If you carry any medical device such as pacemaker and other medical aids on their body or pregnant, do not use the Product.
6. Take off your hair accessories, contact lens and glasses before the use of the Product.
7. If you had an eye surgery such as LASIK, LASEK or lens implant surgery, do not use the Product. If you insist to use it, make sure that you consult with your doctor first.
8. The Product does not substitute a professional's diagnosis or treatment. If necessary, consult with a professional before using the Product.
9. If there is any abnormal symptom such as rash or inflammation, or you had a surgical operation and are in the middle of treatment in the target area, do not use the Product.
10. The heating function has nothing to with hair loss, atopic dermatitis, acne and allergy.
11. The mode setting becomes active once you put on the Product. Otherwise, it could be damaged by exceeding the air pressure limit.
12. If you feel uncomfortable during use, stop the Product immediately and cut off the power. Then, check the Product.
13. For a charging adaptor, use the enclosed DC5V adapter only. Be careful for electrical safety.
14. Turn on the power after putting on the Product. Do not lean on a wall or walk with the Product on.
15. The Product is designed for home use only. Do not use it for commercial or any other (e.g., medical treatment or diagnosis, etc.) purposes.
16. For cleaning, use the sterilized cotton ball with a neutral detergent. Do not EVER try to clean the Product with water.
17. When the heating function is used, if you are dull in feeling temperature, be careful during use.

Problem Solving (1)

Problem Solving

If a problem occurs during use, check the Product as instructed below. Unless it is solved despite the action, cut off the power and contact the Technical Service Center for help.

Technical Service Center : 1588-2336 * Do not disassemble or repair the Product.

Trouble	Cause	Solution
Power is not turned on	Under-charge or discharge	Charge the battery.
Weak air pressure & vibration	Under-charge or discharge	Charge the battery.
Not working	Poor connection	Check if the charger is properly connected.

Problem Solving (2)

Specification

Adapter Voltage	DC 5V
Power Consumption	5W
Dimension / Weight	(W)210mmX(L)254mmX(H)242mm / 1075g
Operating	Temperature : 5℃ ~ 40℃
Conditions	Humidity : ≤80% Atmospheric Pressure : 860hpa ~ 1060hpa

Maintenance

Storage

Keep the Product safe from humidity, rodents and contamination. Before you clean the Product, switch off the power first. Keep the Product in a safe place after use. Do not leave a heavy object on the Product. If a defect is found, contact the engineer. If the Product is left unused for more than 3 months, check if it is normally operated. Then, charge the system and leave it on for about 30 minutes. Do not keep the Product near the wall or roof. Store it in an airy place. Avoid direct sunlight and corrosive gases.

Transportation

The Product is earthquake-resistant to some extent. It is resistant to mild vibrations such as flight, railway, road and ocean shipping. Avoid a strong impact such as snow, rain and electric shock. Make sure that it is positioned in a right direction all the time.

Power Supply

Check the power supply. If the voltage exceeds the rated range, do not turn on the Product. Check the charging unit and cable periodically. If any damage is found, do not use the Product. Send it to the Technical Service Center.

Supplies and Accessories

To get supplies or accessories, visit our website (www.koreabreo.com) or nearby store. If you have any question about them, contact the Customer Service Center at 1588-2336.

Guidance on Quality Assurance

About Warranty

With this warranty, you are able to get the repair & maintenance services easily. Only 1 copy is provided for each product. If lost, it may not be reissued. Be careful not to lose it.

* The warranty period is two years from the date of purchase.

Terms and Conditions

1. You can get repair services free of charge during the warranty period (for two years from the date of purchase). If it is hard to make out the serial number on the Product, or it is not matched with the serial number on your warranty, warranty service may not be available. (A battery is not warranted.)
2. A problem attributable to a user's negligence such as incorrect voltage or fortuitous accident or defect caused by an unqualified repair company will not be warranted.
3. A cost could occur for the service other than warranty services or transportation service.
4. A problem attributable to a user's fault such as incorrect storage or caused by an unqualified engineer or defect which occurred because of force majeure will not be warranted. Furthermore, warranty service is not provided unless the warranty certificate is presented.
5. Enclose the warranty certificate when you apply for the warranty service. Korea Breo is responsible for collecting, processing and recycling the used and disposed products.
6. If used for purposes other than home use (e.g., commercial purpose, etc.), the Product will not be warranted.
7. Any failure caused by natural disasters such as typhoon, earthquake or flood will not be warranted.

* Notice: Please be advised that the services are continuously upgraded so that specific rules are subject to change without a separate notice.

Product Warranty

Product Name :

Serial No. :

Date of Purchase :

Place of Purchase :

Customer Name :

TEL :

Address :

Exchange & Refund Policy

1. If the Product needs be repaired because of a defect arising from a normal use within 10 days from the date of purchase, it can be exchanged or refunded by presenting the warranty or receipt.
2. If the Product needs to be repaired because of a serious problem arising from a normal use within one month from the date of purchase, it can be exchanged or repaired. However, the Product other than the defect should be free of damage. The Product can be exchanged, refunded or repaired provided that the User Manual and other components are kept intact.

For more information about this matter, contact the store from which you bought the Product or our Customer Service Center (1588-2336). We will get back to you soon. A battery is not warranted.

Terms and Conditions

1. A warranty service is available for two years from the date of purchase.
If it is hard to make out the serial number on the Product, or it is not matched with the serial number on your warranty, warranty service may not be available. (A battery is not warranted.)
2. A problem attributable to a user's fault such as the use of incorrect voltage or incorrect storage or caused by an unqualified engineer will not be warranted.
3. A cost could occur for the service other than warranty services or transportation service.
4. A problem attributable to a user's fault such as incorrect storage or caused by an unqualified engineer or defect which occurred because of force majeure will not be warranted. Furthermore, warranty service is not provided unless the warranty certificate is presented.
5. Enclose the warranty certificate when you apply for the warranty service. Korea Breo is responsible for collecting, processing and recycling the used and disposed products.
6. If used for purposes other than home use (e.g., commercial purpose, etc.), the Product will not be warranted.
7. Any failure caused by natural disasters such as typhoon, earthquake or flood will not be warranted.

* Notice: Please be advised that the services are continuously upgraded so that specific rules are subject to change without a separate notice.

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

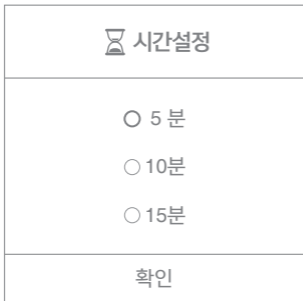
How to Use (1)

1. Please download iBrain relaxer APP
(through Android play store/iOS App Store)
2. Put on iBrain relaxer and start the massage by turning on the power.
3. Open the bluetooth in the mobile phone set and select the "iBrain relaxer" in the bluetooth device list searched.
4. Connect the iBrain relaxer to the bluetooth of the mobile phone.
5. When the product is successfully connected to the bluetooth, the builtin music will stop. (If the built-in music is still in play, please check if the bluetooth of the mobile phone is successfully connected)
6. Run the iBrain relaxer APP installed on the phone.
7. When the connection choosing image appears, connect the iBrain relaxer again in the bluetooth device list searched. If the clues of "can not find the equipment" appears:
 - 1) Please confirm whether the phone bluetooth is connected with the ibrain relaxer
 - 2) Click "find again"

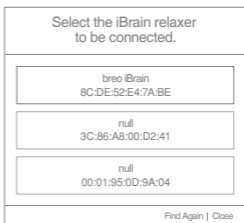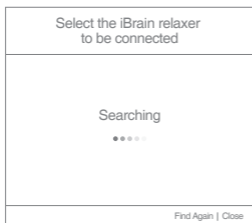8. If the APP bluetooth is not successfully connected, click on the bluetooth button in the APP and select the iBrain relaxer in the bluetooth list appeared.
9. After disconnecting the bluetooth, it takes 5 seconds to reconnect.
10. The iBrain relaxer is successfully connected to the phone bluetooth.
(after connecting the phone bluetooth, the steps of 3-5 could be ignored when reusing the product)
11. Choose the required massage mode.


MESSAGE MODE.	
 mode 1	IF STRESSED →Message your head and neck and ease tense muscles.
 mode 2	TIRED EYES →Get rid of the fatigue in the eyes through massage.
 mode 3	STIFF NECK →Soothe and relax tense muscles in the neck.
 mode 4	HEADACHE →Massage the temple and ease the pain.
 mode 5	I NEED CONCENTRATION →Clear your head through massage.

How to Use (2)


12. Adjust the massage strength.

13. The bluetooth connection and music playback can not be executed under battery symbol flashing state and low electricity condition.
14. Press  and set the timer.

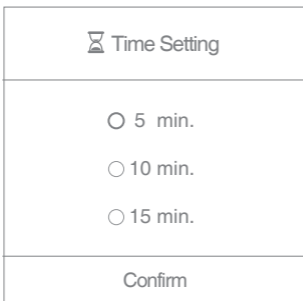15. Press  and set the heating function.

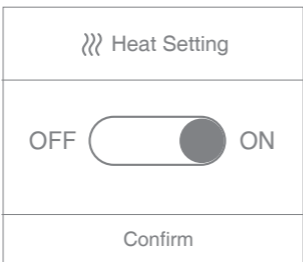16. Press  and set the music.