# K1

User Manual

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# K1

User Manual

Thank you for purchasing  
the 10minds breo products.  
Before using this product,  
please read this manual carefully  
for convenient use.

## Contents

---

Product Overview and Components	27
Product features	28
How to use	29
Specifications and Troubleshooting	30
Maintenance	31
Precautions	32
Replacement and refund	33
Warranty Guide	34
Product Warranty	35

## Product Overview and Components

---

### Product Overview

Manage your knees that bent and stretch all day long with K1 for a lifetime. The ergonomically designed K1 wraps your knees snugly and performs an acupressure treatment on your fine muscles.

### Product Components

K1, User Manual, Charging adapter

K1


User Manual

Charging adapter

## Product Features

### Product Features

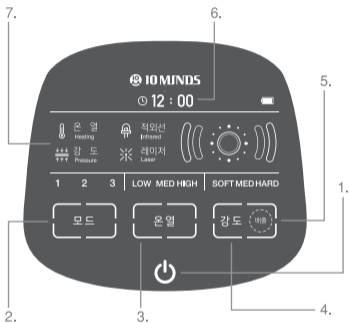1. Power button
2. Mode display
3. Temp. Adjust
4. Pressure Adjust
5. Air pressure discharge
6. Operation Time
7. Functions under currently use

## How to use

---

### Power/Operation

1. Make sure the product is fully charged.
2. Fix the product to your lap and touch the power button to turn on the product.
3. Touch Mode / Heating / Pressure button shortly to use it.
4. The product usage time is automatically set to 15 minutes in total for 1-time use.

### Mode Change

Mode	1	2	3
Heating	O	O	O
Pressure	O	O	O
Infrared (14lights)	O	X	O
Laser (5 lights)	O	O	X

### Mode/Heating/Pressure Change

1. The first mode is set to 1 after power on  
(In order of 1-2-3-2-1)
2. Heating is set to Middle during initial operation.  
(in order of Middle-High-Middle-Low)
3. Pressure is set to Middle during initial operation  
(in order of Middle-High-Middle-Low)

## Specifications and Troubleshooting

### Specifications

Rating Voltage	DC 5V 2.5A	Power Consumption	12.5W
Size/Weight	(W)168mmX(L)193mmX(H)144mm / 900g		
Conditioning Temp	Low 39±3℃, Middle 41±3℃, High 43±3℃		
Use environment	Temp. : 0℃ ~ 40℃ Humid : ≤80% Atmosphere : 860hpa ~ 1060hpa	Laser wavelength	808nm

### Troubleshooting

If you have any problem in using the product, please check the instructions below. If this does not work, turn off the power and send it to the service center. Please note that if you are not an expert, do not disassemble or repair the product personally since we may not provide a further service.

Problem	Possible Cause	Solution
The charger is connected but no light in ON/OFF buttons	Charger is in failure	Replace the charger
LCD touch screen is not activated	Water on the LCD screen	Dry the product surface. Touch the product with dried hand only.

## Maintenance

---

### Storage precautions

Keep the product in a safe place to avoid exposure to moisture and contamination. If the storage period exceeds 6 months, check the operation of product and turn on the power for about 30 minutes. Avoid the direct sunlight and corrosive gases.

### Transporting precautions

This product is designed for seismic resistance and can be transported by plane, train, car, and ship. Avoid snow, rain, and impact as much as possible and put it in the right direction.

### Troubleshooting

If a malfunction occurs during operation, check the main battery first. If the trouble cannot be solved, turn off the power and send it to our A/S service center. Do not self-disassemble or repair.

If the voltage exceeds the rated range of product, do not turn on the product. Check the cables of power supply periodically.

## Precautions

---

- ※ Please read the following cautions carefully and accordingly use the Product.
1. Do not use the product while sleeping.
  2. If you have bleeding in your knees, please do not use.
  3. If you have rheumatoid arthritis, please consult your doctor before use.
  4. If you have trauma, bacterial infection or other specific inflammation, please consult your doctor before use.
  5. Pregnant women, infants and other people with certain skin disorders are not allowed to use the product
  6. Do not use the product for the purpose of diagnosis or treatment of disease.
  7. If you feel any pain, stop operation. Press and hold Pressure/ Discharge button to quickly discharge the air.
  8. Do not let the laser light come in direct contact with your eyes.
  9. When cleaning the unit, use a soft and dry cloth. Make sure the unit is dry before use.
  10. Use only the adapter supplied with the product and pay attention for electrical safety.
  11. Wear the product, turn on the power and start operation.
  12. This product has been developed as a household product, so it should not be used for other commercial uses, treatment or diagnosis.
  13. When cleaning, wipe off the product with the disinfected cotton with neutral detergent. Be careful not to let the product come into contact with water.

## Guidance on Product Replacement / Refund

---

1. If a critical repair due to any performance- or function- related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function- related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have maintained their original conditions before unpacked.

If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center (1588-2336). Batteries are not applicable to any warranty service.

## Guidance on Quality Assurance

---

### Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty service is valid for one year from the date of purchase.

### Quality Assurance Conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty. (Batteries are not applicable to such services.)
2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty. (e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster (e.g. typhoon, earthquake, flood, etc.), no warranty service can be provided.

Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

## Product Warranty

---

Name of Product :

Serial Number :

Date Purchased :

Place Purchased :

Name of Customer :

Contact :

Address :

### Guidance on Product Replacement/Refund

1. If a critical repair due to any performance- or function -related defect under a normal operating environment is required within 10 days from the date purchased, you can request for a product replacement or refund with your purchase receipt.
2. If a critical repair due to any performance- or function -related defect under a normal operating environment is required within 1 month from the date purchased, you can request for a product replacement or repair. However, such request can only be valid if the Product has been damaged under a normal operating environment. In addition, the Product can be replaced, refunded or repaired only if the user manual and all other components have main tained their original conditions before unpacked If you intend to make such request above or have any related inquiry, please contact your buyer or our customer center(1588-2336). Batteries are not applicable to any warranty service.

### Warranty information

With the product warranty, you can easily receive various maintenance services. Only one warranty card is issued for each product. Therefore, you should take extra care not to lose the card since it can be really cumbersome to reissue a new card.

※ The warranty service is valid for one year from the date of purchase.

### Quality Assurance Conditions

1. With the product warranty, you can receive various repair services for one year from the date purchased. However, such services may not be properly provided if the serial number on the Product is not legible, or if such serial number differs from the one registered on the product warranty.  
(Batteries are not applicable to such services.)
2. For any defect due to a user's fault (e.g. use of wrong power, accidental failure, etc.) or a repair service by a non-authorized vendor, no warranty service can be provided.
3. Additional expense may incur for any supplementary services not specified in the product warranty (e.g. delivery service, etc.)
4. For any defect due to a user's fault or negligence, a disassembly service by an unqualified expert, no product warranty, or a natural disaster, no warranty service can be provided.
5. When requesting for any quality assurance service, please enclose your product warranty. We are responsible for collecting, treating and recycling any returned or disposed product.
6. If the Product has been used for any commercial purpose (e.g. business use, public utilities, etc.) other than household uses, no warranty service can be provided.
7. For any defect due to a natural disaster (e.g. typhoon, earthquake, flood, etc.), no warranty service can be provided.

Note: Since quality assurance services get continuously upgraded, detailed regulations may be subject to change without any customer notification.

10minds